



Distinguished enjoyment
Taste the beauty of elegance

THIS MOMENT
I AM HERE ONLY FOR YOU



We are seeking for the 'Glory' in your heart along the way.

Hanging Type



Top Mounted



Cart Type






Bluetooth Music Function:
Relaxing patient's stressful emotion




Detachable spittoon 180°rotatable

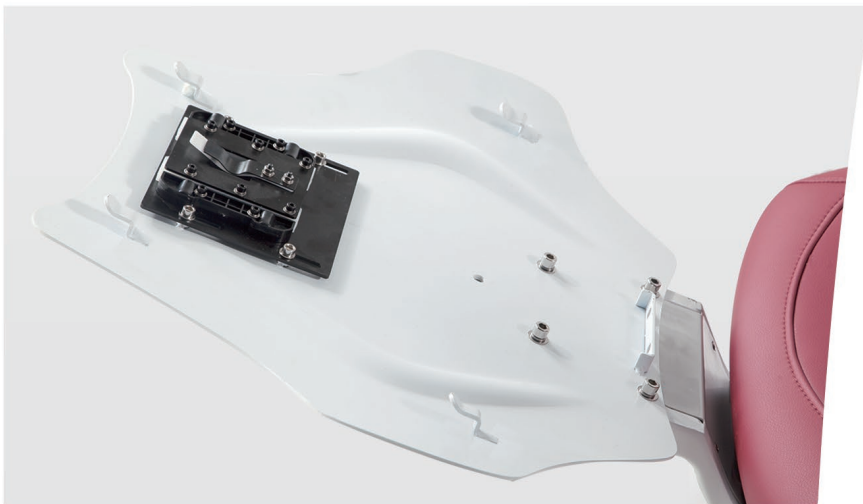

Metal backrest

Detachable backrest convenient for replacement

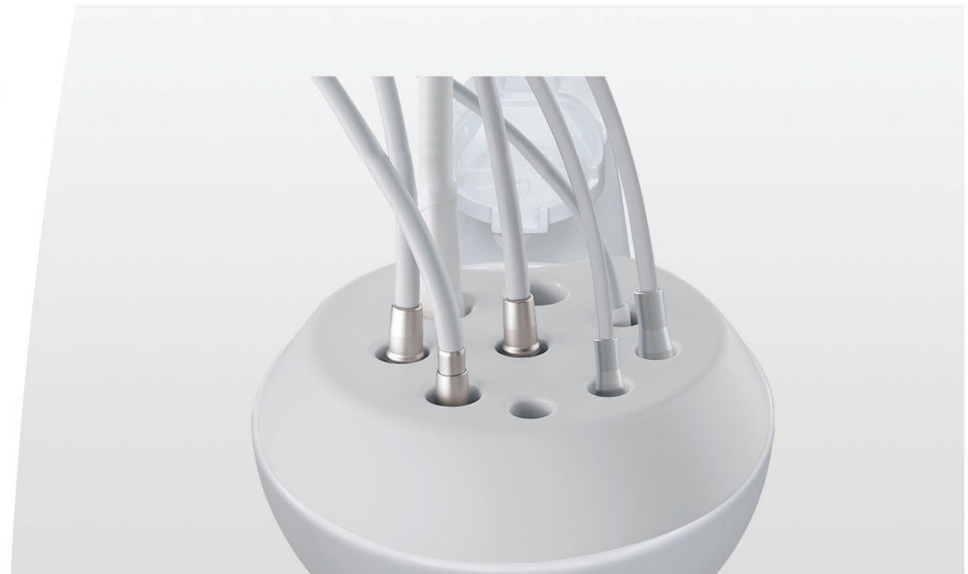

Disinfection

The automatic ASTER system which controlled by one button, provides intensive disinfection for all instrument spray & suction tubing internal circuits with antiseptic liquid




- European style swing type doctor tray.
- Equipped with silicone pad and special tray for instruments.
- Flexible Swing handle, light and easy to operate.

Cingol 3.0 Smart Operating System

- Self- diagram detection system.
- Shock-resistance function maximizes the comfort.
- Clean up chair position, one key to turn off the light, rise the seat to the highest point and the spittoon flushes itself automatically, avoiding blocked pipes and odors, improve work efficiency.
- Safety emergency stop , First-aid chair position.
- Interlock system: when the handpiece or scaler working , the chair movement will lock itself.



- Touch control panel doctor tray.
- Full-ranged acrylic protective cover.
- 135 ° sterilizable silicone handrail cover.
- 3 sets and total 9 memory programs for 3 doctors .
- Spit restoration: the seat and backrest will move to the right position for spitting automatically, and back to the initial position after spitting.
- Water temperature adjustable for all tubes: the water of handpiece, 3-way syringe and scaler tube can be heated.

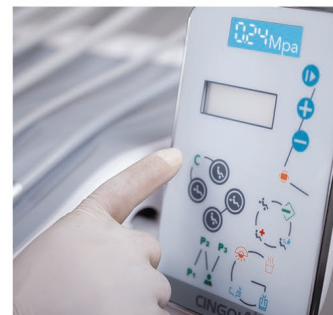

Touch sense control panel



Ti-alloy handpiece tubing connector



Imported-brand 3-way syringe



Multi-joint assistant tray arm.
Multi-angle movement convenient for
four hands operation.



Touch sense assistant tray



3-joint-axes can be positioned for
adults and children or disabled alike

Ergonomics Upholstery

- The specific shape of the ultra-thin backrest design makes more room for the doctor and convenient to operate.
- The flat backrest leaves enough space for doctors & assistance. Legs can freely stretch under the backrest to avoid the soreness caused by V-shaped position during surgery.
- Wear-resistant fibre leather as standard configuration, which is also soft and breathable.
- High elastic sponge ensures the comfort of patients.
- The blockage protection system will stop automatically and rise slightly when the seat encountering resistance during the descent to avoid accidents.
- The backrest and seat cushion move synchronically enables the best movement without back rubbing during the treatment.
- The side box is connected in the middle of the frame, making the whole device more stable.
- Spit restoration: The seat and backrest will move to the right position for spitting automatically, and back to the initial position after spitting.



Sidebox 90°rotatable



▶ Detachable sidebox cover, easy for maintenance.



▶ Rotary plug-in double sealed filter. External design for easy cleaning.




▶ Built-in double bottles water storage system for distilled and disinfected water. The water can pour into the bottles from the top of sidebox, free to detach. Easy to operate and improve the work efficiency.

T.MASTER

European 6-hole
operation light.

Stepless adjustment of shadowless cold light,
close to natural light, Create a more relaxing treatment.

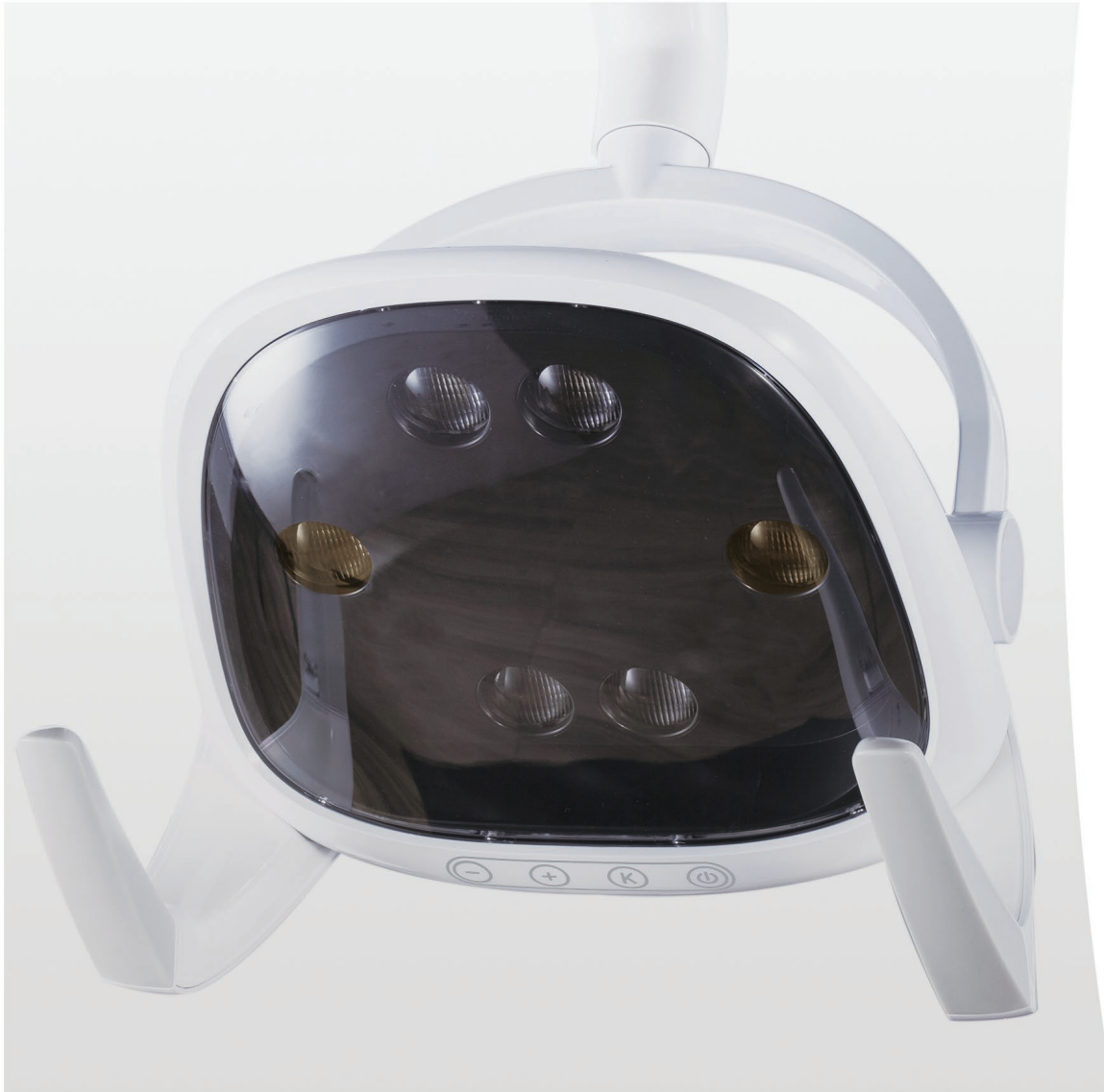

Multi-functional tray
USB charging function



Luxury multi-functional pedal can control the
Operation lamp & Tumbler filling



WINQ brand Italian style hanged type
flexible light arm

OPTIONS



LED curing light



Scaler



Handpiece



Electric motor



Intra oral camera with support



ECG



Binocular loupes & Microscope



T-Master teaching recording lamp



WINQ oilless air compressor



WINQ wet suction




WINQ dry suction




T-Master dental simulator

SPECIFICATION


Input voltage: 230V \pm 10% ,50Hz
 Water pressure:0.2-0.4Mpa/2.0-4.0 bar
 Air pressure: 0.55-0.8Mpa/5.5-8.0 bar
 G.Weight:255Kg

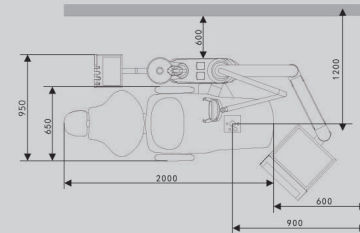

Built-in floor-box installation sketch map

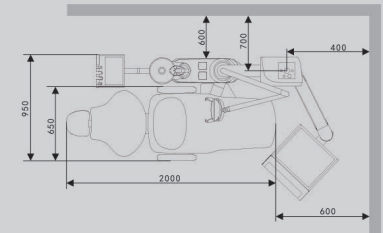

External floor-box installation sketch map

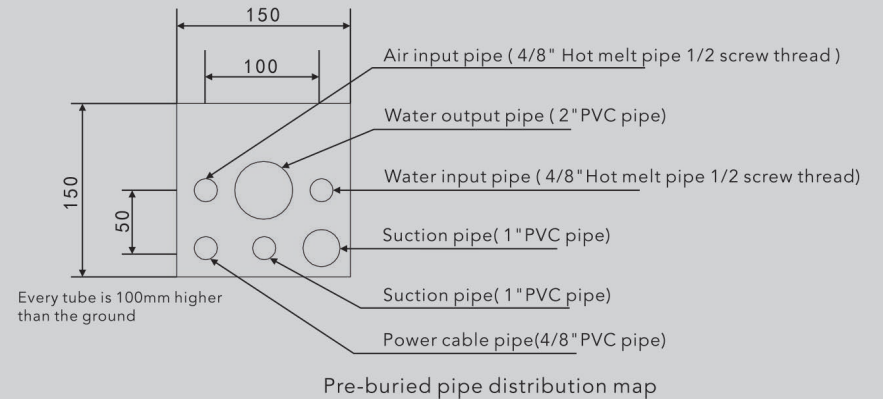

Color Options



> Popular colors ★★★★★

★★★★ C01 ★★★ C03 ★ C10



CINGOL

Concern about International dental Needs,
our Goal is your Obliging Laugh.



Remark:

We reserve the right to improve all the products
without notice in advance;
The product's picture is only for your reference;
the copyright is reserved and belongs to Cingol.

FOSHAN CINGOL MEDICAL INSTRUMENT CO.,LTD
Tel:0086-757-81797099 Fax:0086-757-81797055
Email:cingol@cingol.com QQ:2851731555
Site:www.cingol.com
Add:Building 6-1,Nanshan Industrial Village,Jiangli North Rd,
Danzao Town,Nanhai District,Foshan,China.